# Oral-B®

## Stages POWER

powered by **BRAUN**

Charger Type 3757  
Handle Type 3709

91199111

<b>GB</b>	<b>0 800 731 1792</b>
<b>IE</b>	<b>1 800 509 448</b>
<b>PL</b>	<b>801 127 286</b> <b>801 1 BRAUN</b>
<b>CZ</b>	<b>800 11 33 22</b>
<b>SK</b>	<b>800 333 233</b>
<b>HU</b>	<b>(06-1) 451-1256</b>
<b>HR</b>	<b>01 / 6690 330</b>
<b>SI</b>	<b>080 2822</b>
<b>IL</b>	<b>1-800-666-775</b>



**Internet:**

[www.oralb.com](http://www.oralb.com)


[www.braun.com](http://www.braun.com)

[www.service.braun.com](http://www.service.braun.com)


Braun GmbH  
Frankfurter Straße 145  
61476 Kronberg/Germany/  
Germania/ Германия/ Немачка.


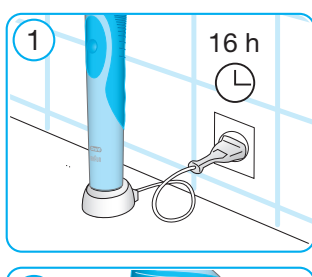
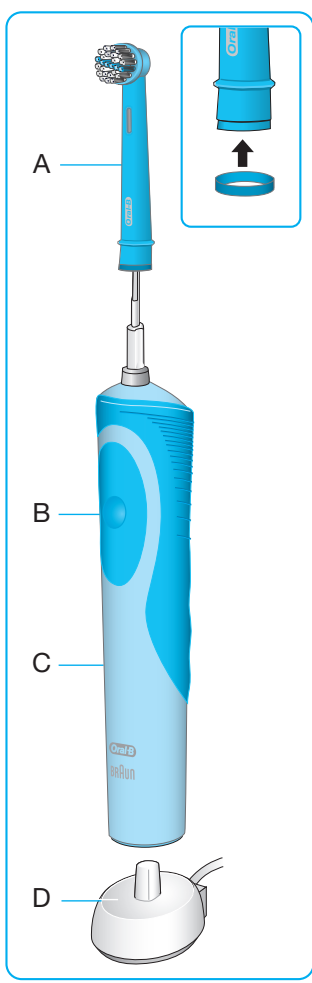
3 outsides

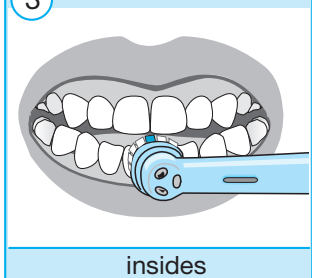insides

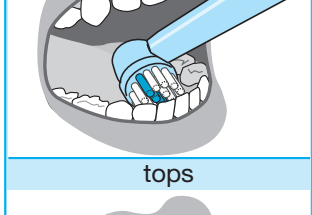tops
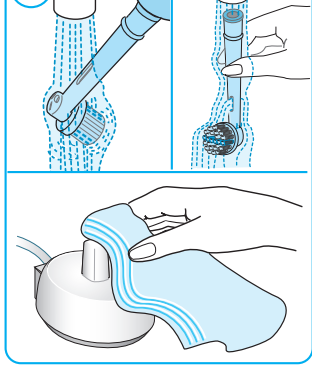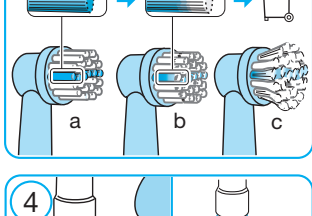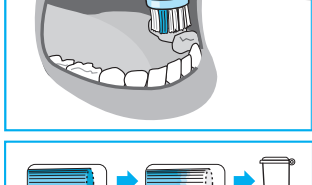
## English

*Dear parents,  
Please read the usage instructions carefully before coaching your child how to use this toothbrush. We recommend that parents supervise their children when they start using the appliance to help them to become familiar with it.*

### **IMPORTANT**

- Periodically check the entire product/cord for damage. A damaged or non-functioning unit should no longer be used. If the product/cord is damaged, take it to an Oral-B Service Centre. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

### **WARNING**

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.

## Description

- A Brush head
- B On/off button
- C Handle
- D Charger

### **Specifications**

For voltage specifications refer to the bottom of the charging unit.

Noise level: ≤65 dB (A)

## Connecting and charging

This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet and place the toothbrush handle (C) on the charging unit. A full charge takes about 16 hours and enables up to five days of regular brushing (twice a day, 2 minutes) (1).
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

## Using the toothbrush

- Wet brush head and apply kids' toothpaste. Guide the brush head to the teeth; then turn on the handle by pressing the on/off switch (B) (2).
- Start with the lower teeth: Slowly guide the brush head from tooth to tooth; first the outsides, then the insides, finally the chewing surfaces (3). Do not press too hard, simply let the brush do all the work.
- After two minutes of operation, an incorporated timer signals by a short stutter the end of the professionally recommended brushing time.
- Turn off the handle while the brush head is still in the mouth.

In the first days of using any electric toothbrush, gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

### Replacement brush heads: Oral-B Stages Power (EB10)

Available at your retailer or at an Oral-B Braun Service Centre.

The brush head is provided with a row of blue INDICATOR® bristles (3a) helping you to monitor brush head replacement need. With proper brushing twice a day for two minutes, the colour will fade halfway in a 3-month period (3b). If bristles begin to splay before the colour starts to recede, this indicates that your child exerts too much pressure on teeth and gums (3c).

## Cleaning recommendations

After brushing, rinse the brush head for several seconds under running water with the handle switched on. Switch off the handle and remove the brush head. Rinse both parts separately under running water then wipe them dry before reassembling the toothbrush. From time to time also clean the charging unit with a damp cloth; never place the charging unit in water (4).

Subject to change without notice.

## Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country

### **Guarantee**

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

## Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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